

GETINGE

GETINGE GS54, GS57
ROTARY SEALERS



A GENERATION OF ROTARY SEALERS: FAST, ENERGY-SAVING AND CONVENIENT.

Getinge's rotary sealers for sterile bags offer lots of attractive features. All models are tilted forward for optimum ergonomics, while at the same time saving space on the packaging table. The GS54 and GS57 sealers have extension options with a extra ergonomic and hygienic packaging table. On both models, the temperature is automatically lowered when the sealer is not in use, which cuts heat emissions and saves energy.

Moreover, all two models GS54 and GS57 achieve a level of performance that no other rotary sealers even come close to. The warm-up time from a cold start is a mere 30 seconds! When the temperature is temporary reduced to save energy, it takes just 20 seconds to regain the working temperature!

Different models for different needs

GS54 is an automatic rotary sealer for smaller clinics. GS57 is fully automatic, with lots of special features. The GS57 also has one or two integrated printers for documentation of process parameters, and can be connected to Getinge's T-DOC® system.

High safety & reliability

Our rotary sealers incorporate a new, thin type of heating element that is highly sensitive to temperature fluctuations. The result is an even temperature and perfect seals.

Tested and CE-marked


Each rotary sealer is tested in accordance with EU safety norms. They are validated and shipped with a certificate, operating instructions, technical information and a list of spare parts.



Easy to change the ink cartridges.



The GS57 can be connected to an external printer (available as an accessory) and computer for printing text on the bag. It can also be linked to a central documentation system.

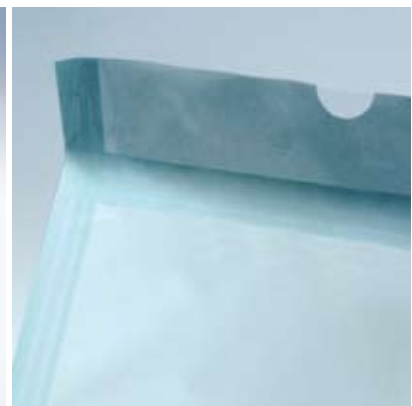

The even temperature assures perfect seals. The GS57 also has one or two integrated printers that can print one or two lines of information on the side of the sterile bag.

GS54



GS57




	GS54	GS57
Width x depth x height, mm	510x280x215	665x280x215
Weight, kg	17	21
Voltage	110-230V/50-60Hz	110-230V/50-60Hz
Wattage	740W	740W
Fuse	4A	4A
Connection	-	USB / Serial

Safety and documentation

GS57/1DK and GS57/2DK allow sealing-pressure control. All GS57 are prepared to be linked with Getinge's T-DOC® system.

Always ready

The heating time from a cold start is approx. 30 seconds. After a pause it's just 20 seconds!

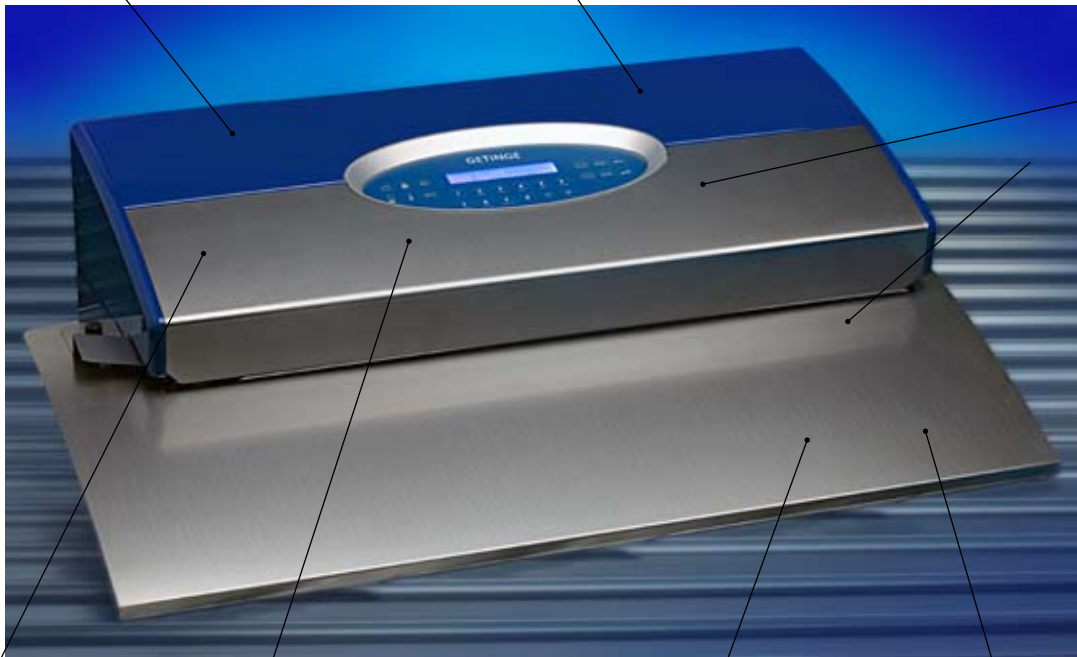

Easy-to-change Ink cartridges

Attractive design

Low heat emission

Our GS54 and GS57 models automatically reduce the temperature when the machine is not being used.

Packing table

As a option a packing table can be added to both models, a big advantage is the hygienic stainless steel plate.

Space saving

Both models requires very little space on the packing table without the optional packing table.



Getinge provides complete solutions for effective and efficient cleaning, disinfection and sterilization in the healthcare and life science sectors. Our know-how comprises everything from architectural planning, production and handling equipment, to systems for full traceability of sterile goods. Our commitment covers expert advice, training and long-term technical support.

GETINGE

Getinge Infection Control
PO Box 69, SE-310 44 Getinge, Sweden
Phone: +46 35 15 55 00
Fax: +46 35 549 52
info@getinge.com www.getinge.com

GETINGE

THE GETINGE GROUP is a leading global provider of equipment and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. Equipment, services and technologies are supplied under the brands **ARJO** for patient hygiene, patient handling and wound care, **GETINGE** for infection control and prevention within healthcare and life science and **MAQUET** for surgical workplaces, cardiopulmonary and critical care.