

GETINGE

GETINGE 1200 FRONT-LOADED FLUSHER-DISINFECTOR





FLUSHER-DISINFECTORS THAT FEEL RIGHT AT HOME

The Getinge 1200 flusher-disinfector proves that hi-tech disinfection equipment doesn't have to look institutional. The gentle contours of the soft blue polymer front panels fit right into every environment where disinfection is required. Yet there is absolutely no compromise on performance.

The outstanding performance of Getinge's stainless steel flusher-disinfectors is well known in the healthcare sector all over the world.

Getinge 1200 is a further step in the development - a flusher-disinfector that combines good design with superior performance in preventing the spread of infection.

Small & easy to position

The Getinge 1200 is compact and just 45 cm wide, so it fits in easily, even where space is limited.

If vertical space is limited, there is also a convenient under-the-counter model.



The Getinge 1200 is available in two different heights: 870 mm (under-the-counter model) and 1255 mm.

SUPERIOR TECHNOLOGY FOR CLEANING & DISINFECTION

When it comes to effective cleaning and disinfection, Getinge goes further than anyone else. In the Getinge 1200, heavily soiled items are first rinsed with water via no fewer than nine fixed and two rotating jet nozzles.


The rinsed items are then steam-disinfected via the same nozzles, which spray directly into every single item. This means that internal disinfection of bottles, etc. is always assured.




The superior performance of the Getinge 1200 not only disinfects the items in the machine. The machine itself -- chamber, pipes, nozzles, etc. -- is also disinfected in every cycle, thanks to PSD (Pipe System Disinfection).

Easy to use

The Getinge 1200 is extremely easy to use for healthcare staff. Just place the item in the holder, close the door and the item is emptied automatically. Then press the start button. After 3-7 minutes the item will be clean, disinfected and ready for use once more!

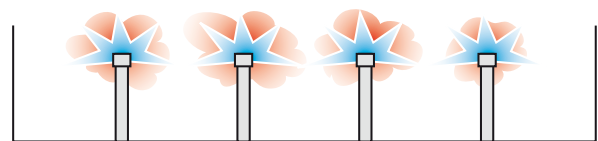

In the Getinge 1200, every nozzle supplies steam, thus injecting steam directly into every item to be cleaned.

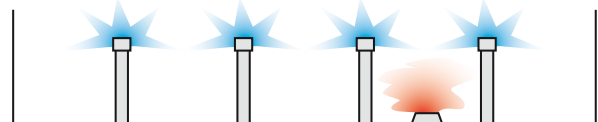

In many other flusher-disinfectors, only one of the nozzles supplies steam, which then has to make its own way into bottles, etc.



User-friendly control panel

The display shows

- current temperature, throughout the process
- fault code in the event of a breakdown

Program for bedpans and toilet bowls

Program for urine bottles

Easy to service

Service is easy to perform from the front, thanks to the extendible electronic box.

Environment-friendly

Satisfies strict environmental requirements: low water consumption, good heat insulation, low noise level and recyclable material.

Safe disinfection

Steam is sprayed in through all nozzles, assuring disinfection of all surfaces of the items, as well as of the machine chamber, pipes and nozzles.

Complies with European standard EN 15883

Hygienic

Smooth edges on the inside of the door and a slender chamber with rounded corners minimize the risk of residual dirt and bacterial growth.

Safe & complete process

The door is closed throughout the entire process to assure staff safety, as well as a thorough cleaning and disinfection process.

Safe storage

For detergents and rinsing agent.





SUPERIOR CLEANING AND DISINFECTION



Suitable for most items

The Getinge 1200 is ideal for cleaning and disinfecting most of the items used in geriatric care.

Simply place the item in the holder (smaller bowls in a basket), close the door, and the item is emptied automatically.

Nozzles for superior cleaning

In the Getinge 1200, both water and steam are supplied via nine fixed nozzles and two rotating ones (one at the top and one in the rear chamber wall).

The result is superior cleaning and disinfection.



Getinge provides complete solutions for effective and efficient cleaning, disinfection and sterilization in the healthcare and life science sectors. Our know-how comprises everything from architectural planning, production and handling equipment, to systems for full traceability of sterile goods. Our commitment covers expert advice, training and long-term technical support.

GETINGE

Getinge Infection Control
PO Box 69, SE-310 44 Getinge, Sweden
Phone: +46 35 15 55 00
Fax: +46 35 18 14 50
info@getinge.com www.getinge.com

GETINGE

THE GETINGE GROUP is a leading global provider of equipment and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. Equipment, services and technologies are supplied under the brands **ARJO** for patient hygiene, patient handling and wound care, **GETINGE** for infection control and prevention within healthcare and life science and **MAQUET** for surgical workplaces, cardiopulmonary and critical care.