

# GETINGE

## GETINGE 1345 WORKTOP WASHER-DISINFECTOR



# GETINGE 1345 – FROM USE TO RE-USE IN 40 MINUTES!

Getinge 1345 is a fully automatic countertop washer-disinfector that minimizes the time between use, cleaning, disinfection and re-use of goods. It's small, so it's easy to place close to contaminated goods. It's available with a built-in drying function. And the overall processing time is short.

Perfect for all small wards and clinics, e.g. dental and district clinics, that need to quickly re-use soiled instruments and circulation goods.



The small footprint enables placing the Getinge 1345 on a standard worktop (60 cm deep). This gives you the perfect working height, as well as easy installation.

## **Built-in drying**

Thanks to the integrated drying function, no manual effort is needed to dry the goods. This also reduces the risk of recontamination. To keep the total process time as short as possible, the drying cycle starts the second the washing is completed. To further speed up drying, the lid also opens slightly.

The drying phase is stopped if the lid is opened fully.

## **Short process time**


Getinge 1345 can be directly connected to the hot water supply and 3-phase electrical current. The total

process time is then just 40 minutes, including drying.

In the standard version, the machine is connected to a cold water supply and 1-phase electrical current, in which case the total process time is 60 minutes.

## **Even for tubular instruments**

For washing the insides of tubular instruments, a flush ramp can be mounted against/towards the back of the chamber.

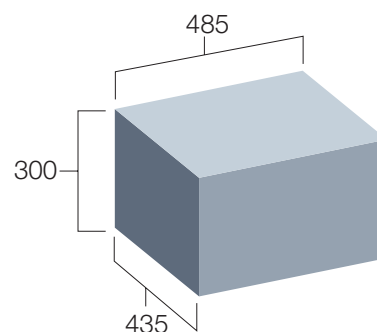



**Fast processes**

Instruments are ready for re-use within 40/60 minutes, including drying.

**Easy placement**

Small size (HxWxD): 505 x 500 x 643 mm.

**Built-in drying (optional)**

No manual labor required.

**High loading capacity**

Large chamber: 300 x 435 x 485 mm.

**Ergonomical**

Can be placed at a convenient working height.

**Printer (optional)**

For process documentation.

**Independent temperature monitoring (optional)**

For maximum process reliability (in accordance with EN ISO 15883).

Different loading alternatives with different goods.





Getinge provides complete solutions for effective and efficient cleaning, disinfection and sterilization in the healthcare and life science sectors. Our know-how comprises everything from architectural planning, production and handling equipment, to systems for full traceability of sterile goods. Our commitment covers expert advice, training and long-term technical support.

# GETINGE

Getinge Infection Control  
PO Box 69, SE-310 44 Getinge, Sweden  
Phone: +46 35 15 55 00  
Fax: +46 35 549 52  
info@getinge.com www.getinge.com

## GETINGE

**THE GETINGE GROUP** is a leading global provider of equipment and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. Equipment, services and technologies are supplied under the brands **ARJO** for patient hygiene, patient handling and wound care, **GETINGE** for infection control and prevention within healthcare and life science and **MAQUET** for surgical workplaces, cardiopulmonary and critical care.